

The ePOSTURO – is a modern computer based posturography system designed for daily work. Center of gravity (COG) is an important component of balance and should be assessed when evaluating someone’s posture. With ePOSTURO software and the eVNG /eHIT mask or the mobile sensor Motion BLUE, static and dynamic Posturography can easily be done. The software calculates all relevant parameters automatically.




BioMed Jena GmbH
 Am Egelsee 1
 D-07743 Jena
 Germany
 Phone: +49-3641 - 35690 -0
 Fax: +49-3641 - 35690 -9
 email: info@biomed-jena.de
 Internet: www.biomed-jena.de

Static and dynamic posturography - ePOSTURO



- » Bluetooth 2.1
- » Wireless charging
- » also with eHIT / eVNG mask
- » +/- 0.5° vertical and 2° horizontal
- » 100 Hz sampling rate
- » animated 3D

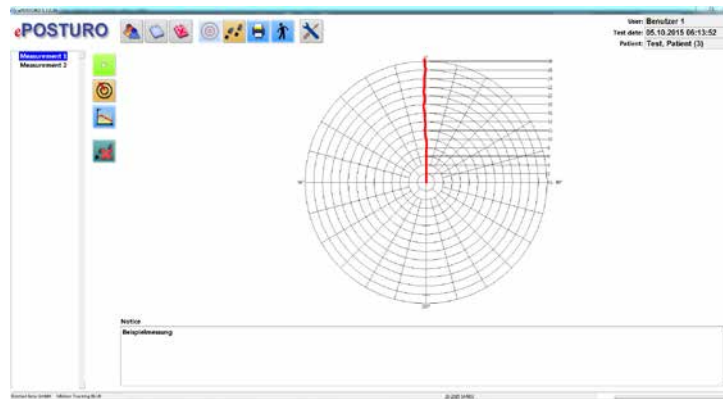

ePOSTURO

Posturography system

- » networking
- » easy to use
- » automatic Unterberger stepping test
- » Risk to fall level by Prof. Walther
- » Treatment suggestions by Prof. Walther
- » innovative result presentation
- » Software: Windows 10/11 with intel Core i5 Full HD recommended
- » 650 mAh battery





Die Posturographie (Gleichgewichtsanalyse) ist ein Verfahren zur Ermittlung der Funktionsfähigkeit der Gleichgewichtsregulation. Mit dem neuartigen ePOSTURO ist es möglich, die Bewegung des Körpers im Raum in Echtzeit zu messen.

Die Software ermittelt automatisch alle interessanten Parameter.



BioMed Jena GmbH
Am Egelsee 1
D-07743 Jena
Germany
Phone: +49-3641 - 35690 -0
Fax: +49-3641 - 35690 -9
email: info@biomed-jena.de
Internet: www.biomed-jena.de




Statische & dynamische Posturographie ePOSTURO

ePOSTURO

Posturographie-System

- » Bluetooth 2.1
- » kontaktloses Laden
- » eHIT/eVNG Maske als Sensor
- » max. 3 synchron aufzeichnende Sensoren
- » +/- 0.5° vertikal and 2° horizontal
- » 100 Werte / s
- » animierte 3D Darstellung

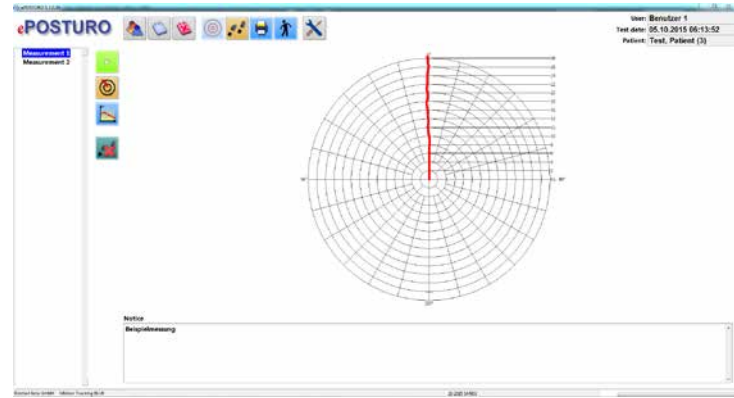

- » Netzwerkfähig
- » Easy to use
- » automatischer Unterberger-Tretversuch
- » Sturzrisiko-Bewertung nach Walther
- » Trainingsempfehlungen nach Walther
- » Romberg-Quotient
- » Kompatibel mit Windows 10/11

